

Drying and heating chambers



100

Avantgarde.Line ED | FD | FP | FED series

Classic.Line E | ED | FD | FED | M series

Solid.Line ED-S | FD-S series



Drying and heating chambers Avantgarde.Line

Series ED Avantgarde.Line	Page 6
Series FD Avantgarde.Line	Page 12
Series FED Avantgarde.Line	Page 19
Series FP Avantgarde.Line	Page 26

BINDER drying chambers with natural convection or forced convection are renowned for their high quality and reliability. Thanks to their extensive redevelopment, these all-rounders are now even more firmly focused on addressing the requirements involved in drying or heating samples. What's more, in addition to the best energy efficiency on the market and convenient operation, the latest APT.line[™] technology offers outstanding temperature accuracy to give the Avantgarde.Line incubators real appeal.

PRODUCT SELECTION DRYING AND HEATING CHAMBERS AVANTGARDE. LINE

> Page 6	› Page 12	Series FED Avantgarde. Line > Page 19	Series FP Avantgarde.Line > Page 26
•	•	•	٠
•	•	٠	•
•	•	•	•
•	•	•	•
300	-	_	-
_	300	300	300
_	300	300	300
	• • • 300	- 300	- - - - 300 - - - 300 300 300

 $\bullet Standard \ equipment$

ADVANTAGES AT A GLANCE



EASY DATA LOGGING

Series ED, FD, FED

The USB connector, standard on the ED, FD and FED, provides easy data logging. Ethernet connectivity is standard on FED and an option for ED and FD.



BEST ENERGY EFFICIENCY ON THE MARKET AND LOW INSTALLATION REQUIREMENTS

Series ED, FD, FED


thanks to low heat dissipation and an optimized heating system.



EASY, INTUITIVE HANDLING

Series ED, FD, FED

Our state-of-the-art controller with LCD display, excellent readability, alphanumeric display, and various functions provide a user friendly interface. Choice of languages as standard: English, German, and French.

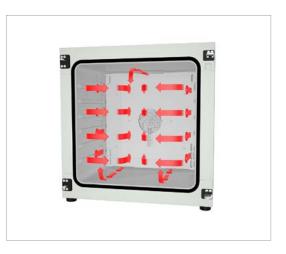

COMPLETELY HOMOGENEOUS TEMPERATURE DISTRIBUTION

Series FD, FED

This feature provides excellent spatial and temporal temperature accuracy with the all-new generation of APT.line™ technology.

Series ED Avantgarde.Line | with natural convection

The strengths of a BINDER ED series drying chamber include routine drying and sterilization tasks up to 300 °C. Thanks to the natural convection, all thermal processes which use this drying chamber are highly efficient. The ED series ensures fast, uniform drying.



ED 56 model

BENEFITS

- · Best energy efficiency on the market
- Excellent temporal and spatial temperature accuracy
- USB connection for recording data
- Temperature range up to +300 °C



MAIN FEATURES

- Temperature range: +5 °C above ambient temperature to +300 °C
- Up to 30% lower energy consumption compared to conventional units on the market
- High temperature accuracy thanks to APT.line[™] technology
- Natural convection
- Controller with LCD display
- Electromechanical control of the exhaust air flap
- 2 chrome-plated racks
- Units up to 115 liters are stackable
- Class 2 integrated independent adjustable temperature safety device (DIN 12880) with visual alarm
- Ergonomic handle design

AVAILABLE SIZES

Model	Interior volume [L]	Online data sheets
ED 56	57	› go2binder.com/en-ED56-AL
ED 115	114	› go2binder.com/en-ED115-AL
ED 260	255	› go2binder.com/en-ED260-AL
ED 720	743	› go2binder.com/en-ED720-AL

G

OPTIONAL EQUIPMENT

Access port with silicone plug	For introducing external measuring instruments. Large selection of different diameters and positions (left, right, or top)	
Viewing window and interior lighting	Viewing window in the door with a triple-glazed insulating glass and interior lighting for easy visibility of the samples.	
Shelf racks	Perforated shelves or shelf racks, chrome-plated or made from stainless steel.	
HEPA air filter	Particle filter on fresh air supply for unit; filter classification H14 (acc. to. EN 1822-1:2009, > 99.995% at 0.3 μm).	
Services	You can access our installation, maintenance, calibration, validation, or warranty services by taking out a Maintenance Service Agreement or by simply getting in touch when you need us. More information is available in the accessories and services chapter.	
BINDER INDIVIDUAL	Using a BINDER series production product as a basis, we develop individual solutions as a single unit or in small batches. All products are tested, certified, and supplied with a full functional guarantee and warranty. More information is available in the accessories and services chapter.	- 22
► All extras online go2bind	er.com/en-options	

Series ED Avantgarde.Line | Drying and heating chambers Avantgarde.Line

TECHNICAL DATA


Description	ED 56	ED 115	ED 260	ED 720
Measures				
Interior volume [L]	57	114	255	743
Net weight of the unit (empty) [kg]	42	57	85	169
Permitted load [kg]	70	150	270	315
Load per rack [kg]	30	30	40	45
Wall clearance back [mm]	160	160	160	160
Wall clearance sidewise [mm]	100	100	100	100
Internal Dimensions				
Width [mm]	360	510	610	960
Height [mm]	420	530	760	1280
Depth [mm]	380	425	550	605
Housing dimensions not incl. fittings and connections				
Width net [mm]	560	710	810	1165
Height net [mm]	625	705	940	1590
Depth net [mm]	565	605	760	870
Performance Data Temperature				
Temperature range +5 °C above ambient temperature to [°C]	300	300	300	300
Temperature variation at 150 °C [± K]	2.3	1.5	1.8	3.2
Temperature fluctuation at 150 °C [± K]	0.4	0.4	0.8	0.8
Heating-up time to 150 °C [min]	50	45	50	85
Recovery time after 30 seconds door open at 150 °C [min]	19	16	20	25
Tabellen				
Air change (approx.) at 100 °C [x/h]	8	6	3	2
Electrical data				
Nominal power (version 120 V) [kW]	1.15	1.35	_	-
Nominal power (version 230 V) [kW]	1.05	1.25	2.25	-
Nominal power (version 240 V) [kW]	-	-	2.45	-
Nominal power (version 400 V) [kW]	-	-	_	4.1
Fixtures				
Number of shelves (std./max.)	2/4	2/5	2/8	2/16

ORDERING INFORMATION [Art. No.]


Further details: Go to www.binder-world.com > Search > Enter Article Number

Option model	ED 56	ED 115	ED 260	ED 720
Standard	-	_	_	9010-0341
Standard	-	-	9010-0340	-
Standard	9010-0333	9010-0335	9010-0339	-
Standard	9010-0334	9010-0336	-	-
	Standard Standard Standard	Standard – Standard – Standard 9010-0333	Standard - - Standard - - Standard - - Standard 9010-0333 9010-0335	Standard - - Standard - - Standard - 9010-0333 9010-0335 Standard 9010-0333 9010-0335 9010-0339

DIMENSIONS incl. fittings and connections [mm]

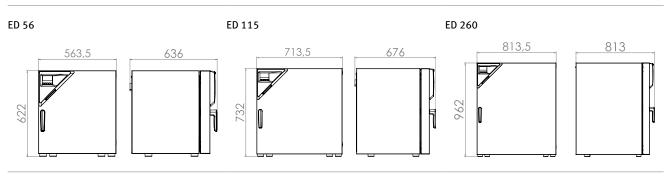

ED 720

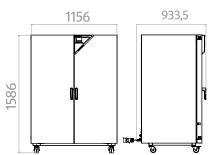

ORDERING INFORMATION OPTIONS

Designation	Description	ED 56	ED 115	ED 260	ED 720	Art. No.
	back					
	10 mm	•	•	•	•	8012-128
	30 mm	•	•	•	•	8012-054
	50 mm	•	•	•	•	8012-054
	left					
	10 mm	•	•	•	•	8012-127
	30 mm	•	•	•	•	8012-101
	50 mm	•	•	•	•	8012-102
	100 mm	-	•	•	•	8012-103
ccess port with sili- one plug	right					
	10 mm	•	•	•	•	8012-126
	30 mm	•	•	•	•	8012-101
	50 mm	•	•	•	•	8012-102
	100 mm	_	•	•	•	8012-103
	top					
	10 mm	•	•	•	•	8012-126
	30 mm	•	•	•	•	8012-053
	50 mm	•	•	•	•	8012-054
	100 mm	•	•	•	•	8012-102
Alarm function for overheating	Switchable acoustic alarm, with adjustable limit value on the indepen- dent temperature safety device	•	•	•	•	8012-163
Analog output 4-20 mA	for temperature values (output not adjustable)	•	•	•	•	8012-162
Calibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test temperature	•	•	•	•	8012-111

Designation	Description	ED 56	ED 115	ED 260	ED 720	Art. No.
	for temperature, measurement in center of chamber at specified tempe- rature	•	•	•	•	8012-113
Calibration certificate,	temperature measurement incl. certificate and 27 measuring points at specified temperature	•	•	•	•	8012-159
temperature	temperature measurement incl. certificate, 15- 18 measuring points at specified temperature		•	•	•	8012-157
	temperature measurement incl. certificate, 9 measuring points at speci- fied temperature	•	•	•	•	8012-155
Class 3.1 independent temperature safety device	with visual alarm (DIN 12880)	•	•	•	•	8012-098
Clock	battery-backed	•	•	•	•	8012-164
		•	-	-	-	8012-100
D	made of FKM, silicone-free		•	-	-	8012-100
Door gasket			-	•	-	8012-101
		-	-	-	•	8012-101
		•	•	•	-	8012-166
Door lock	lockable door handle		-	-	•	8012-181
Tabaan at intanfaaa	for Multi Management Coffiger ADT COMTM	•	•	•	-	8012-164
Ethernet interface	for Multi Management Software APT-COM™	-	-	-	•	8012-183
Inner chamber, rein- forced	max. total load 350 kg	-	-	-	•	8012-182
	with 15 W light bulb					
	120 V option model	•	•	-	-	8012-163
Interior lighting	230/240 V option model	•	•	•	-	8012-162
	with two 15 W light bulbs					
	400 V option model	-	-	-	•	8012-182
Measurement of air exchange rate	in accordance with ASTM D5374, definition and protocol according to ambient temperature	•	•	•	•	8012-019
Pt 100 temperature sensor	additional flexible Pt 100, interior, for displaying the temperature on the unit display	•	•	•	•	8012-161
	viewing window (triple insulating glass) in door, 240 x 240 mm, and 15 W interior					
	120 V option model	•	_	-	-	8012-105
	230 V option model	•	-	-	-	8012-105
	viewing window (triple insulating glass) in door, 300 x 480 mm, and 30 W interior					
Viewing window and	400 V option model	-	-	-	•	8012-183
interior lighting	viewing window (triple insulating glass) in door, 370 x 370 mm, and 15 W interior					
	120 V option model	-	•	-	-	8012-106
		-	•	-	-	8012-106
	230 V option model	-	-	•	-	8012-106
	240 V option model	_	_	•	_	8012-106

ORDERING INFORMATION ACCESSORIES

Designation	Description	ED 56	ED 115	ED 260	ED 720	Art. No.
APT-COM™ 4 BA- SIC-Edition	for simple logging and documentation requirements with up to 5 networ- ked units.					
SIC-Euron	version 4, BASIC edition	•	•	•	•	9053-0039

Designation	Description	ED 56	ED 115	ED 260	ED 720	Art. No.
APT-COM™ 4 GLP-Editi- on	for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11.					
	version 4, GLP edition	•	•	•	•	9053-004
APT-COM TM 4 PROFES-	convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units.					
SIONAL-Edition	version 4, PROFESSIONAL edition	•	•	•	•	9053-004
Ethernet interface	for Multi Management Software APT-COM™	•	•	•	-	8012-164
oH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	•	•	•	•	8012-225
	IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature, CO_2 , O_2 , pressure, depending on unit					
	Digital in PDF format	•	•	•	•	7057-0003
Qualification docu-	Hard copy inside folder	•	•	•	•	7007-000
nents	IQ/OQ/PQ documents – supporting documents for validation performed by customers, according to customer requirements, PQ section added to qualification folder IQ/OQ; parameters: temperature, CO ₂ , O ₂ – or pres- sure, depending on unit					
	Digital in PDF format	•	•	•	•	7057-000
	Hard copy inside folder	•	•	•	•	7007-000
			-	-	-	8012-203
	chrome plated	-	•	-	-	8012-204
		-	-	•	-	8012-204
		_	-	-	•	8012-204
Rack		•	-	-	_	8012-217
		_	•	-	-	8012-2172
	stainless steel	-	-	•	-	8012-217
		_	-	-	•	8012-2174
		_	-	•	-	8012-218
Rack, heavy load	Stainless steel, max. load per rack 70 kg	-	-	-	•	8012-218
Rubber pads	set anti-slip feet	•	•	•	_	8012-203
		•	-	-	-	8012-217
		-	•	-	-	8012-218
Shelf, perforated	Stainless steel	-	-	•	-	8012-218
		-	-	-	•	8012-218
	stable cart, casters with locking brakes, dimensions: W 1,000 x D 800 x H 780 mm	•	•	_	_	9051-001
Table on castors	stable cart, casters with locking brakes, dimensions: W 1,300 x D 800 x H 780 mm	_	_	•	_	9051-001

Series FD Avantgarde.Line | with forced convection

A BINDER FD series Avantgarde.Line heating oven is always used when fast drying and sterilization is required. Thanks to its fully homogeneous temperature distribution, quick dynamics and powerful fan, this heating oven saves valuable time.



FD 56 model

BENEFITS

- Best energy efficiency on the market
- Excellent temporal and spatial temperature accuracy
- USB connection for recording data
- Temperature range up to +300 °C



MAIN FEATURES

- Temperature range: +10 °C above ambient temperature to +300 °C
- High temperature accuracy thanks to APT.line[™] technology
- Forced convection
- Controller with LCD display
- Electromechanical control of the exhaust air flap
- 2 chrome-plated racks
- Units up to 115 L are stackable
- Class 2 independent temperature safety device (DIN 12880) with visual temperature alarm
- Ergonomic handle design
- USB port for recording data

AVAILABLE SIZES

Model	Interior volume [L]	Online data sheets
FD 56	60	› go2binder.com/en-FD56-AL
FD 115	116	› go2binder.com/en-FD115-AL
FD 260	259	> go2binder.com/en-FD260-AL
FD 720	741	> go2binder.com/en-FD720-AL

G

OPTIONAL EQUIPMENT

Access port with silicone plug	For introducing external measuring instruments. Large selection of different diameters and positions (left, right, or top)	
Viewing window and interior lighting	Viewing window in the door with a triple-glazed insulating glass and interior lighting for easy visibility of the samples.	1
Shelf racks	Perforated shelves or shelf racks, chrome-plated or made from stainless steel.	
HEPA air filter	Particle filter on fresh air supply for unit; filter classification H14 (acc. to. EN 1822-1:2009, > 99.995% at 0.3 μm).	
Services	You can access our installation, maintenance, calibration, validation, or warranty services by taking out a Maintenance Service Agreement or by simply getting in touch when you need us. More information is available in the accessories and services chapter.	
BINDER INDIVIDUAL	Using a BINDER series production product as a basis, we develop individual solutions as a single unit or in small batches. All products are tested, certified, and supplied with a full functional guarantee and warranty. More information is available in the accessories and services chapter.	-

All extras online go2binder.com/en-options

Series FD Avantgarde.Line | Drying and heating chambers Avantgarde.Line

TECHNICAL DATA


Description	FD 56	FD 115	FD 260	FD 720
Measures				
Interior volume [L]	60	116	259	741
Net weight of the unit (empty) [kg]	41	57	85	166
Permitted load [kg]	70	150	270	315
Load per rack [kg]	30	30	40	45
Wall clearance back [mm]	160	160	160	160
Wall clearance sidewise [mm]	100	100	100	100
Internal Dimensions				
Width [mm]	400	550	650	1000
Height [mm]	440	550	780	1300
Depth [mm]	345	385	515	570
Housing dimensions not incl. fittings and connections				
Width net [mm]	560	710	810	1165
Height net [mm]	625	735	940	1590
Depth net [mm]	565	605	825	870
Performance Data Temperature				
Temperature range +10 °C above ambient temperature to [°C]	300	300	300	-
Temperature range +12 °C above ambient temperature to [°C]	_	-	-	300
Temperature variation at 150 °C [± K]	1.7	1.7	1.9	2.5
Temperature fluctuation at 150 °C [± K]	0.3	0.3	0.4	0.6
Heating-up time to 150 °C [min]	15	18	19	25
Recovery time after 30 seconds door open at 150 °C [min]	4	4	5	6
Tabellen				
Air change (approx.) at 100 °C [x/h]	80	32	9	5
Electrical data				
Nominal power (version 120 V) [kW]	1.2	1.4	-	-
Nominal power (version 230 V) [kW]	1.1	1.3	2.3	-
Nominal power (version 240 V) [kW]	-	-	2.5	-
Nominal power (version 400 V) [kW]	-	-	-	4.5
Fixtures				
Number of shelves (std./max.)	2/4	2/5	2/8	2/16

ORDERING INFORMATION [Art. No.]


Further details: Go to www.binder-world.com > Search > Enter Article Number

Voltage	Option model	FD 56	FD 115	FD 260	FD 720
400 V 3~ 50/60 Hz	Standard	-	_	-	9010-0311
240 V 1~ 60 Hz	Standard	-	-	9010-0310	-
230 V 1~ 50/60 Hz	Standard	9010-0303	9010-0305	9010-0309	-
120 V 1~ 60 Hz	Standard	9010-0304	9010-0306	-	-

DIMENSIONS incl. fittings and connections [mm]

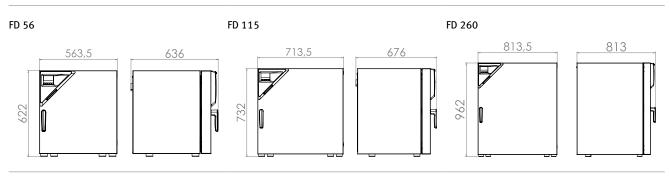

FD 720

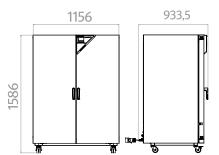

ORDERING INFORMATION OPTIONS

Designation	Description	FD 56	FD 115	FD 260	FD 720	Art. No.
	left					
	10 mm	•	•	•	•	8012-127
	30 mm	•	•	•	•	8012-101
	50 mm	•	•	•	•	8012-102
	100 mm	-	•	•	•	8012-102
	right					
	10 mm	•	•	•	•	8012-126
Access port with sili- cone plug	30 mm	•	•	•	•	8012-101
	50 mm	•	•	•	•	8012-101
	100 mm	-	•	•	•	8012-102
	top					
	10 mm	•	•	•	•	8012-126
	30 mm	•	•	•	•	8012-136
	50 mm	•	•	•	•	8012-138
	100 mm	•	•	•	•	8012-167
Alarm function for overheating	Switchable acoustic alarm, with adjustable limit value on the indepen- dent temperature safety device	•	•	•	•	8012-163
Analog output 4-20 mA	for temperature values (output not adjustable)	•	•	•	•	8012-162
Calibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test temperature	•	•	•	•	8012-111

Designation	Description	FD 56	FD 115	FD 260	FD 720	Art. No.
	for temperature, measurement in center of chamber at specified tempe- rature	•	•	•	•	8012-113
Calibration certificate,	temperature measurement incl. certificate and 27 measuring points at specified temperature	•	•	•	•	8012-159
emperature	temperature measurement incl. certificate, 15- 18 measuring points at specified temperature	•	•	•	•	8012-157
	temperature measurement incl. certificate, 9 measuring points at speci- fied temperature	•	•	•	•	8012-155
Class 3.1 independent emperature safety levice	with visual alarm (DIN 12880)		•	•	•	8012-164
Clock	battery-backed	•	•	•	•	8012-164
		•	-	-	-	8012-164
		_	•	-	_	8012-165
)oor gasket	made of FKM, silicone-free	_	-	•	_	8012-165
		_	_	_	•	8012-183
		•	•	•	_	8012-166
oor lock	lockable door handle	_	_	_	•	8012-181
		•	•	•	_	8012-09
thernet interface	for Multi Management Software APT-COM™	_	_	_	•	8012-09
an, reinforced	for increased air change rates, incl. uprated heat output for faster heating up time and improved recovery time; max. 300 °C; not in con- junction with load management option					
	400 V option model	-	-	_	•	8012-10
		•	_	-	_	8012-168
as-tight option model	Incl. additional measures for the greatest possible impermeability to gas; not in conjunction with access port or interior lighting option	-	•	-	-	8012-16
		-	-	•	_	8012-16
	for faster heating up time and improved recovery time; max. 300 °C					
eating output, rein-	400 V option model	_	-	-	•	8012-10
prced	120/230 V option model	•	_	_	_	8012-10
		_	•	_	_	8012-10
		•	_	_	_	8012-10
	On fresh air supply for unit; filter classification H14 (acc. to. EN 1822-	_	•	_	_	8012-10
EPA air filter	1:2009, > 99.995% at 0.3 µm); not in conjunction with reinforced heating	_	_	•	_	8012-10
	output, reinforced fan		_	_	•	8012-10
		•	_	_	_	8012-16
ert gas connection	With gas inlet and outlet, Ø 10 mm, incl. additional measures for the greatest possible impermeability to gas; not in conjunction with access	_	•	_	_	8012-16
left gas connection	port or interior lighting option		_	•	_	8012-16
	may total load 170 kg		•	_		
	max. total load 170 kg			•		8012-10
nner chamber, rein- orced	max. total load 300 kg	-	-	-	-	8012-10
, iccu	max. total load 350 kg	-	-	-	•	8012-18
	max. total load 85 kg	•	-	-	-	8012-10
	with 15 W light bulb					
	120 V option model	•	•	-	-	8012-16
Interior lighting	230/240 V option model	•	•	•	-	8012-162
	with two 15 W light bulbs					
	400 V option model	-	-	-	•	8012-182
leasurement of air xchange rate	in accordance with ASTM D5374, definition and protocol according to ambient temperature	•	•	•	•	8012-119

Designation	Description	FD 56	FD 115	FD 260	FD 720	Art. No.
	viewing window (triple insulating glass) in door, 240 x 240 mm, and 15 W interior					
	120 V option model	•	-	-	-	8012-170
	230 V option model	•	_	-	-	8012-170
	viewing window (triple insulating glass) in door, 300 x 480 mm, and 30 W interior					
Viewing window and	400 V option model	-	_	-	•	8012-183
interior lighting	viewing window (triple insulating glass) in door, 370 x 370 mm, and 15 W interior					
	120 V option model	-	•	-	-	8012-170
	220 Mantian madel	-	•	-	-	8012-170
	230 V option model	-	_	•	-	8012-170
	240 V option model	_	-	•	_	8012-171

ORDERING INFORMATION ACCESSORIES

Designation	Description	FD 56	FD 115	FD 260	FD 720	Art. No.
APT-COM™ 4 BA- SIC-Edition	for simple logging and documentation requirements with up to 5 networ- ked units.					
SIC-Edition	version 4, BASIC edition	•	•	•	•	9053-003
APT-COM™ 4 GLP-Editi- on	for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11.					
	version 4, GLP edition	•	•	•	•	9053-004
APT-COM™ 4 PROFES-	convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units.					
SIONAL-Edition	version 4, PROFESSIONAL edition	•	•	•	•	9053-004
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	•	•	•	•	8012-2250
	IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature, CO ₂ , O ₂ , pressure, depending on unit					
	Digital in PDF format	•	•	•	•	7057-000
Qualification docu-	Hard copy inside folder	•	•	٠	•	7007-000
ments	IQ/OQ/PQ documents – supporting documents for validation performed by customers, according to customer requirements, PQ section added to qualification folder IQ/OQ; parameters: temperature, CO_2 , O_2 – or pres- sure, depending on unit					
	Digital in PDF format	•	•	•	•	7057-000
	Hard copy inside folder	•	•	•	•	7007-000
		•	-	-	-	8012-2038
		-	•	-	-	8012-204
	chrome plated	-	-	•	-	8012-2042
		-	-	-	•	8012-204
Rack		•	-	-	-	8012-208
		-	•	-	-	8012-216
	stainless steel	-	_	•	-	8012-2169
		-	-	-	•	8012-2170
		-	-	•	-	8012-218
Rack, heavy load	Stainless steel, max. load per rack 70 kg	-	-	-	•	8012-218
Rubber pads	set anti-slip feet	•	•	•	_	8012-203

Series FD Avantgarde.Line | Drying and heating chambers Avantgarde.Line

Designation	Description		FD 115	FD 260	FD 720	Art. No.
Shelf, perforated		•	-	-	-	8012-2175
		-	•	-	-	8012-2176
	Stainless steel	-	-	•	-	8012-2177
		-	-	-	•	8012-2178
Table on castors	stable cart, casters with locking brakes, dimensions: W 1,000 x D 800 x H 780 mm	•	•	-	-	9051-0018
	stable cart, casters with locking brakes, dimensions: W 1,300 x D 800 x H 780 mm	_	_	•	_	9051-0019

Series FED Avantgarde.Line | with forced convection and enhanced timer functions

BINDER heating chambers of the FED Avantgarde.Line series provide almost unlimited capacity and are particularly adaptable. With advanced timing functions and a controllable fan, temperature and convection conditions are easily controlled.



FED 56 model

BENEFITS

- Best energy efficiency on the market
- Excellent temporal and spatial temperature accuracy
- USB connection for recording data
- Temperature range up to 300 °C



MAIN FEATURES

- Temperature range: +10 °C above ambient temperature to +300 °C
- Up to 30% lower energy consumption compared to conventional units on the market
- High temperature accuracy thanks to APT.line[™] technology
- Adjustable fan speed
- Controller with LCD display and enhanced time functions
- Electromechanical control of the exhaust air flap
- 2 chrome-plated racks
- Units up to 115 liters are stackable
- Class 2 integrated independent adjustable temperature safety device (DIN 12880) with visual alarm
- Ergonomic handle design

AVAILABLE SIZES

G

Model	Interior volume [L]	Online data sheets
FED 56	60	› go2binder.com/en-FED56-AL
FED 115	116	› go2binder.com/en-FED115-AL
FED 260	259	› go2binder.com/en-FED260-AL
FED 720	741	› go2binder.com/en-FED720-AL

OPTIONAL EQUIPMENT

Access port with silicone plug	For introducing external measuring instruments. Large selection of different diameters and positions (left, right, or top)	
Viewing window and interior lighting	Viewing window in the door with a triple-glazed insulating glass and interior lighting for easy visibility of the samples.	1
Shelf racks	Perforated shelves or shelf racks, chrome-plated or made from stainless steel.	
HEPA air filter	Particle filter on fresh air supply for unit; filter classification H14 (acc. to. EN 1822-1:2009, \rightarrow 99.995% at 0.3 μm).	
Services	You can access our installation, maintenance, calibration, validation, or warranty services by taking out a Maintenance Service Agreement or by simply getting in touch when you need us. More information is available in the accessories and services chapter.	
BINDER INDIVIDUAL	Using a BINDER series production product as a basis, we develop individual solutions as a single unit or in small batches. All products are tested, certified, and supplied with a full functional guarantee and warranty. More information is available in the accessories and services chapter.	- 22


► All extras online go2binder.com/en-options


TECHNICAL DATA

Description	FED 56	FED 115	FED 260	FED 720
Measures				
Interior volume [L]	60	116	259	741
Net weight of the unit (empty) [kg]	41	57	85	166
Permitted load [kg]	70	150	270	315
Load per rack [kg]	30	30	40	45
Wall clearance back [mm]	160	160	160	160
Wall clearance sidewise [mm]	100	100	100	100
Internal Dimensions				
Width [mm]	400	550	650	1000
Height [mm]	440	550	780	1300
Depth [mm]	345	385	515	570
Housing dimensions not incl. fittings and connections				
Width net [mm]	560	710	810	1165
Height net [mm]	625	705	940	1590
Depth net [mm]	565	605	760	870
Performance Data Temperature				
Temperature range +10 °C above ambient temperature to [°C]	300	300	300	-
Temperature range +12 °C above ambient temperature to [°C]	-	-	-	300
Temperature variation at 150 °C [± K]	1.4	1.2	1.6	2.0
Temperature fluctuation at 150 °C [± K]	0.3	0.3	0.4	0.6
Heating-up time to 150 °C [min]	15	18	19	25
Recovery time after 30 seconds door open at 150 °C [min]	4	4	5	6
Electrical data				
Nominal power (version 120 V) [kW]	1.2	1.4	-	-
Nominal power (version 208 V) [kW]	-	-	-	4.5
Nominal power (version 230 V) [kW]	1.1	1.3	2.3	-
Nominal power (version 240 V) [kW]	-	-	2.5	-
Nominal power (version 400 V) [kW]	-	-	-	4.5
Fixtures				
Number of shelves (std./max.)	2/4	2/5	2/8	2/16


ORDERING INFORMATION [Art. No.] Further details: Go to www.binder-world.com > Search > Enter Article Number

Voltage	Option model	FED 56	FED 115	FED 260	FED 720
400 V 3~ 50/60 Hz	Standard	-	-	_	9010-0301
240 V 1~ 60 Hz	Standard	-	-	9010-0300	-
230 V 1~ 50/60 Hz	Standard	9010-0295	9010-0293	9010-0299	-
208 V 3~ 60 Hz	Standard	-	-	-	9010-0302
120 V 1~ 60 Hz	Standard	9010-0296	9010-0294	-	-

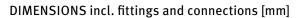

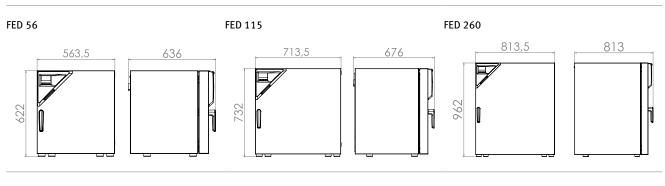

FED 720

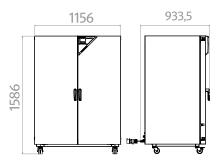

ORDERING INFORMATION OPTIONS

Designation	Description	FED 56	FED 115	FED 260	FED 720	Art. No.
	left					
	10 mm	•	•	•	•	8012-127
	30 mm	•	•	•	•	8012-165
	50 mm	•	•	•	•	8012-165
	100 mm	-	•	•	•	8012-166
	right					
	10 mm	•	•	•	•	8012-127
Access port with sili- cone plug	30 mm	•	•	•	•	8012-165
one prus	50 mm	•	•	•	•	8012-165
	100 mm	-	•	•	•	8012-167
	top					
	10 mm	•	•	•	•	8012-126
	30 mm	•	•	•	•	8012-136
	50 mm	•	•	•	•	8012-138
	100 mm	•	•	•	•	8012-167
Alarm function for overheating	Switchable acoustic alarm, with adjustable limit value on the indepen- dent temperature safety device	•	•	•	•	8012-163
Analog output 4-20 mA	for temperature values (output not adjustable)	•	•	•	•	8012-162
Calibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test temperature	•	•	•	•	8012-111

Designation	Description	FED 56	FED 115	FED 260	FED 720	Art. No.
	for temperature, measurement in center of chamber at specified tempe- rature	•	•	•	•	8012-1137
Calibration certificate,	temperature measurement incl. certificate and 27 measuring points at specified temperature	•	•	•	•	8012-1596
temperature	temperature measurement incl. certificate, 15- 18 measuring points at specified temperature	•	•	•	•	8012-157
	temperature measurement incl. certificate, 9 measuring points at speci- fied temperature	•	•	•	•	8012-1554
Class 3.1 independent temperature safety device	with visual alarm (DIN 12880)	•	•	•	•	8012-164
Clock	battery-backed	•	•	•	•	8012-164
		•	_	-	_	8012-164
- • ·		-	•	-	-	8012-165
Door gasket	made of FKM, silicone-free	_	-	•	-	8012-165
		_	-	-	•	8012-183
		•	•	•	_	8012-166
Door lock	lockable door handle	_	_	-	•	8012-181
Fan, reinforced	for increased air change rates, incl. uprated heat output for faster heating up time and improved recovery time; max. 300 °C; not in con- junction with load management option					
	400 V option model	-	-	-	•	8012-178
Gas-tight option model		•	_	-	-	8012-168
	Incl. additional measures for the greatest possible impermeability to gas; not in conjunction with access port or interior lighting option	-	•	-	-	8012-168
		-	-	•	-	8012-168
	for faster heating up time and improved recovery time; max. 300 °C					
Heating output, rein- forced	120/230 V option model	•	_	-	_	8012-169
		-	•	-	-	8012-169
	On fresh air supply for unit; filter classification H14 (acc. to. EN 1822- 1:2009, > 99.995% at 0.3 μm); not in conjunction with reinforced heating output, reinforced fan	•	-	-	-	8012-169
		_	•	-	-	8012-169
HEPA air filter		_	_	•	_	8012-169
		_	_	-	•	8012-183
		•	_	_	_	8012-169
Inert gas connection	With gas inlet and outlet, Ø 10 mm, incl. additional measures for the greatest possible impermeability to gas; not in conjunction with access	_	•	_	_	8012-169
C	port or interior lighting option	_	_	•	_	8012-169
	max. total load 170 kg	_	•	_	_	8012-171
Inner chamber, rein-	max. total load 300 kg	_	_	•	_	8012-171
forced	max. total load 350 kg	_	_	_	•	8012-182
	max. total load 85 kg	•	_	_	_	8012-171
	with 15 W light bulb					
	120 V option model	•	•	_	_	8012-163
	230/240 V option model	•	•	•	_	8012-162
Interior lighting	with two 15 W light bulbs	-				
	208 V option model	_	_	_	•	8012-097
	400 V option model	_	_	_	•	8012-182
Measurement of air exchange rate	in accordance with ASTM D5374, definition and protocol according to ambient temperature	•	•	•	•	8012-182
Pt 100 temperature sensor	additional flexible Pt 100, interior, for displaying the temperature on the unit display	•	•	•	•	8012-162

Series FED Avantgarde.Line | Drying and heating chambers Avantgarde.Line

Designation	Description	FED 56	FED 115	FED 260	FED 720	Art. No.
	viewing window (triple insulating glass) in door, 240 x 240 mm, and 15 W interior					
	120 V option model	•	-	-	-	8012-1703
	230 V option model	•	-	-	-	8012-1701
	viewing window (triple insulating glass) in door, 300 x 480 mm, and 30 W interior					
Viewing window and	208/400 V option model	-	-	-	•	8012-1832
nterior lighting	viewing window (triple insulating glass) in door, 370 x 370 mm, and 15 W interior					
	120 V option model	-	•	-	-	8012-1702
		-	•	-	_	8012-170
	230 V option model	-	-	•	-	8012-170
	240 V option model	_	_	•	_	8012-171

ORDERING INFORMATION ACCESSORIES

Designation	Description	FED 56	FED 115	FED 260	FED 720	Art. No.
APT-COM™ 4 BA- SIC-Edition	for simple logging and documentation requirements with up to 5 networ- ked units.					
SIC-Edition	version 4, BASIC edition	•	•	•	•	9053-0039
APT-COM™ 4 GLP-Editi- on	for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11.					
	version 4, GLP edition	•	•	•	•	9053-0042
APT-COM™ 4 PROFES-	convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units.					
SIONAL-Edition	version 4, PROFESSIONAL edition	•	•	•	•	9053-0040
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	•	•	•	•	8012-2250
	IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature, CO ₂ , O ₂ , pressure, depending on unit					
Qualification docu- ments	Digital in PDF format	•	•	•	•	7057-0001
	Hard copy inside folder	•	•	•	•	7007-0001
	IQ/OQ/PQ documents – supporting documents for validation performed by customers, according to customer requirements, PQ section added to qualification folder IQ/OQ; parameters: temperature, CO ₂ , O ₂ – or pres- sure, depending on unit					
	Digital in PDF format	•	•	•	•	7057-0005
	Hard copy inside folder	•	•	•	•	7007-0005
		•	-	-	-	8012-2038
		_	•	-	-	8012-2040
	chrome plated	_	-	•	-	8012-2042
		_	-	-	•	8012-2044
Rack		•	-	-	-	8012-2085
		-	•	-	-	8012-2168
	stainless steel	_	-	•	_	8012-2169
		_	-	-	•	8012-2170
		-	-	-	•	8012-2170 8012-2184
Rack, heavy load	Stainless steel, max. load per rack 70 kg	-				

Designation	Description	FED 56	FED 115	FED 260	FED 720	Art. No.
				200	,20	
		•	-	-	-	8012-2175
Shelf, perforated	Stainless steel	-	•	-	-	8012-2176
		-	-	•	-	8012-2177
		-	-	-	•	8012-2178
Table on castors	stable cart, casters with locking brakes, dimensions: W 1,000 x D 800 x H 780 mm	•	•	-	-	9051-0018
	stable cart, casters with locking brakes, dimensions: W 1,300 x D 800 x H 780 mm	-	-	•	-	9051-0019

Series FP Avantgarde.Line | with forced convection and program functions

A BINDER material test chamber with mechanical convection of the FP series provides reliably short drying times and particularly fast heating – even for chambers under full loads.





FP 115 model

BENEFITS

- Uniform drying conditions thanks to APT.line[™] technology
- Identical test conditions throughout the chamber interior independent of sample size and quantity
- Outstanding thermal insulation saves operation costs
- Advanced programming options

MAIN FEATURES

- Temperature range: +10 °C above ambient temperature up to +300 °C
- High temperature accuracy thanks to APT.line[™] technology
- Adjustable fan speed
- Controller with interval and real-time programming
- Electromechanical control of the exhaust air flap
- 2 chrome-plated racks
- Units up to 115 liters are stackable
- Class 2 independent adjustable temperature safety device (DIN 12880) and class 3.1 independent temperature safety device (DIN 12880) with visual/acoustic alarm available for selection on the controller
- Ergonomic handle design
- Computer interface: Ethernet

AVAILABLE SIZES

Model	Interior volume [L]	Online data sheets
FP 56	60	› go2binder.com/en-FP56
FP 115	116	› go2binder.com/en-FP115
FP 260	259	> go2binder.com/en-FP260
FP 720	741	> go2binder.com/en-FP720

G

OPTIONAL EQUIPMENT

Access port with silicone plug	For introducing external measuring instruments. Large selection of different diameters and positions (left, right, or top)	
Viewing window and interior lighting	Viewing window in the door with a triple-glazed insulating glass and interior lighting for easy visibility of the samples.	
Shelf racks	Perforated shelves or shelf racks, chrome-plated or made from stainless steel.	
HEPA air filter	Particle filter on fresh air supply for unit; filter classification H14 (acc. to. EN 1822-1:2009, > 99.995% at 0.3 μm).	
Services	You can access our installation, maintenance, calibration, validation, or warranty services by taking out a Maintenance Service Agreement or by simply getting in touch when you need us. More information is available in the accessories and services chapter.	
BINDER INDIVIDUAL	Using a BINDER series production product as a basis, we develop individual solutions as a single unit or in small batches. All products are tested, certified, and supplied with a full functional guarantee and warranty. More information is available in the accessories and services chapter.	-

► All extras online go2binder.com/en-options

Series FP Avantgarde.Line | Drying and heating chambers Avantgarde.Line

TECHNICAL DATA

Description	FP 56	FP 115	FP 260	FP 720
Measures				
Interior volume [L]	60	116	259	741
Net weight of the unit (empty) [kg]	44	62	96	194
Permitted load [kg]	70	150	270	315
Load per rack [kg]	30	30	40	45
Wall clearance back [mm]	160	160	160	160
Wall clearance sidewise [mm]	100	100	100	100
Internal Dimensions				
Width [mm]	400	550	650	1000
Height [mm]	400	550	780	1300
Depth [mm]	342	383	512	574
Housing dimensions not incl. fittings and connections				
Width net [mm]	560	710	810	1160
Height net [mm]	622	732	962	1590
Depth net [mm]	596	636	786	816
Performance Data Temperature				
Temperature range +10 °C above ambient temperature to [°C]	300	300	300	-
Temperature range +12 °C above ambient temperature to [°C]	-	-	-	300
Temperature variation at 150 °C [± K]	1.4	1.4	1.7	2.4
Temperature fluctuation at 150 °C [± K]	0.4	0.3	0.4	0.6
Heating-up time to 150 °C [min]	15	19	19	25
Recovery time after 30 seconds door open at 150 °C [min]	4	5	6	6
Tabellen				
Air change (approx.) at 100 °C [x/h]	80	32	9	5
Electrical data				
Nominal power (version 120 V) [kW]	1.2	1.4	2.5	_
Nominal power (version 208 V) [kW]	-	-	-	4.5
Nominal power (version 230 V) [kW]	1.2	1.4	2.5	-
Nominal power (version 400 V) [kW]	-	-	-	4.5
Fixtures				
Number of shelves (std./max.)	2/5	2/5	2/8	2/16

ORDERING INFORMATION [Art. No.]

Further details: Go to www.binder-world.com > Search > Enter Article Number

Voltage	Option model	FP 56	FP 115	FP 260	FP 720
400 V 3~ 50/60 Hz	Standard	-	_	-	9010-0366
230 V 1~ 50/60 Hz	Standard	9010-0363	9010-0364	9010-0365	-
208 V 3~ 60 Hz	Standard	-	-	-	9010-0370
120 V 1~ 60 Hz	Standard	9010-0367	9010-0368	9010-0369	_

Designation	Description	FP 56	FP 115	FP 260	FP 720	Art. No.
	left					
	10 mm	•	•	•	•	8012-219
	30 mm	•	•	•	•	8012-221
	50 mm	•	•	•	•	8012-221
	100 mm	_	•	•	•	8012-220
	right					
	10 mm	•	•	•	•	8012-220
Access port with sili-	30 mm	•	•	•	•	8012-221
cone plug	50 mm	•	•	•	•	8012-22
	100 mm	_	•	•	•	8012-22
	top					
	10 mm	•	•	•	•	8012-21
	30 mm	•	•	•	•	8012-22
	50 mm	•	•	•	•	8012-22
	100 mm	•	•	•	•	8012-22
Analog output 4-20 mA	for temperature values (output not adjustable)	•	•	•	•	8012-21
Calibration certificate,	for temperature; for extending the measurement in center of chamber to include another test temperature	•	•	•	•	8012-22
Calibration certificate, temperature	for temperature, measurement in center of chamber at specified tempe- rature	•	•	•	•	8012-22
	temperature measurement incl. certificate and 27 measuring points at specified temperature	•	•	•	•	8012-22
	temperature measurement incl. certificate, 15- 18 measuring points at specified temperature	•	•	•	•	8012-22
	temperature measurement incl. certificate, 9 measuring points at speci- fied temperature	•	•	•	•	8012-22
	made of FKM, silicone-free	•	-	-	-	8012-22
		-	•	-	-	8012-22
Door gasket		-	-	•	-	8012-22
		-	-	-	•	8012-22
		•	•	•	-	8012-22
Door lock	lockable door handle	_	_	-	•	8012-22
		•	_	-	_	8012-22
Gas-tight option model	Incl. additional measures for the greatest possible impermeability to	_	•	-	_	8012-22
	gas; not in conjunction with access port or interior lighting option	_	_	•	_	8012-22
		•	_	_	_	8012-22
	On fresh air supply for unit; filter classification H14 (acc. to. EN 1822-	_	•	_	_	8012-22
IEPA air filter	1:2009, > 99.995% at 0.3 μm); not in conjunction with reinforced heating output, reinforced fan	_	_	•	_	8012-22
	output, remoteed fan	_	_	_	•	8012-22
		•	_	_	_	8012-21
nert gas connection	With gas inlet and outlet, Ø 10 mm, incl. additional measures for the greatest possible impermeability to gas; not in conjunction with access	_	•	_		8012-21
	port or interior lighting option	_	_	•	_	8012-21
	max. total load 170 kg	_	•	_	_	8012-21
				•		
nner chamber, rein- orced	max. total load 300 kg	-	-	-	•	8012-22
	max. total load 350 kg	•				8012-22
Measurement of air	max. total load 85 kg in accordance with ASTM D5374, definition and protocol according to	-	-	-	•	8012-22

ORDERING INFORMATION OPTIONS

Designation	Description	FP 56	FP 115	FP 260	FP 720	Art. No.
Pt 100 temperature sensor	additional flexible Pt 100, interior, for displaying the temperature on the unit display	•	•	•	•	8012-2193

ORDERING INFORMATION ACCESSORIES

Designation	Description	FP 56	FP 115	FP 260	FP 720	Art. No.
APT-COM™ 4 BA-	for simple logging and documentation requirements with up to 5 networ- ked units.					
SIC-Edition	version 4, BASIC edition	•	•	•	•	9053-0039
APT-COM™ 4 GLP-Editi- on	for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11.					
	version 4, GLP edition	•	•	•	•	9053-004
APT-COM™ 4 PROFES- SIONAL-Edition	convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units.					
SIONAL-Edition	version 4, PROFESSIONAL edition	•	•	•	•	9053-004
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	•	•	•	•	8012-225
	IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature, CO_2 , O_2 , pressure, depending on unit					
	Digital in PDF format	•	•	•	•	7057-000
Qualification docu-	Hard copy inside folder	•	•	•	•	7007-000
ments	IQ/OQ/PQ documents – supporting documents for validation performed by customers, according to customer requirements, PQ section added to qualification folder IQ/OQ; parameters: temperature, CO_2 , O_2 – or pres- sure, depending on unit					
	Digital in PDF format	•	•	•	•	7057-000
	Hard copy inside folder	•	•	•	•	7007-000
		•	-	-	-	8012-203
		-	•	-	-	8012-204
	chrome plated	_	-	•	-	8012-204
		-	-	-	•	8012-204
Rack		•	-	-	-	8012-208
		-	•	-	-	8012-216
	stainless steel	-	-	•	-	8012-216
		-	-	-	•	8012-217
Rack accessories	fasteners (1 set of 4) for additional security of racks	•	•	•	•	8012-053
		_	-	•	-	8012-218
Rack, heavy load	Stainless steel, max. load per rack 70 kg	-	-	-	•	8012-218
Rubber pads	set anti-slip feet	•	•	•	-	8012-203
		•	_	_	_	8012-217
		-	•	-	-	8012-217
Shelf, perforated	Stainless steel	-	-	•	-	8012-217
		_	_	-	•	8012-217
	stable cart, casters with locking brakes, dimensions: W 1,000 x D 800 x H 780 mm	•	•	_	_	9051-001
Table on castors	stable cart, casters with locking brakes, dimensions: W 1,300 x D 800 x H 780 mm	_	_	•	_	9051-001

Drying and heating chambers Classic.Line

Series E Classic.Line	Page 36
Series ED Classic.Line	Page 40
Series FD Classic.Line	Page 45
Series FED Classic.Line	Page 50
Series M Classic.Line	Page 55

With natural or forced convection, all BINDER drying and heating chambers offer first-rate quality and process reliability. With its wide temperature range from 5 °C above ambient temperature to 300 °C, the heating chamber can be used in all applications. Short heating-up times and large power reserves are ideal for applications in science, research, and industry. The excellent thermal insulation also reduces operating costs.


In its FP and M series, BINDER offers material test chambers with forced convection and customized programming.

PRODUCT SELECTION DRYING AND HEATING CHAMBERS CLASSIC.LINE

Features	Series E Classic.Line > Page 36	Series ED Classic. Line > Page 40	Series FD Classic. Line > Page 45	Series FED Classic. Line Page 50	Series M Classic.Line > Page 55
Model					
Model 23	-	•	٠	-	_
Model 28	•	_	_	-	_
Model 53	_	_	_	_	•
Model 115	_	_	_	_	•
Model 240	-	-	_	-	•
Model 400	_	•	_	•	•
Model 720	_	-	_	-	•
General data					
Temperature range [°C]	60230	-	_	-	_
Temperature range +5 °C above ambient temperature to [°C]	_	300	300	300	300

Standard equipment

ADVANTAGES AT A GLANCE

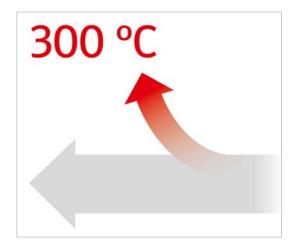

BROAD RANGE OF APPLICATIONS

Series ED, FD, FED, M

Many different applications can take advantage of the wide temperature range from 5 °C above ambient temperature to 300 °C. The FP and M series also impress with their adjustable air change rates and large power reserves.



HOMOGENEOUS TEMPERATURE CONDITIONS WITH NATURAL CONVECTION

Series ED


Gentle temperature distribution provided by the APT.line™ preheating chamber technology.



HOMOGENEOUS TEMPERATURE CONDITIONS WITH FORCED CONVECTION

Series FD, FED, M

Long-term operational stability even under full load from the APT.line[™] preheating chamber technology.

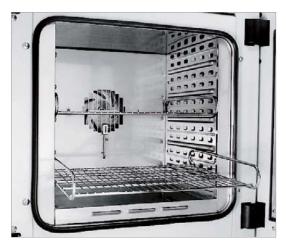

CONVENIENT WORK ENVIRONMENT

Series E, FD, FED, M

BINDER interior concept: highlights rack with tilt protection, no permanent fixtures and stain-less-steel inner chamber.

Series E Classic.Line | with mechanical adjustment

With a BINDER E series heating oven, you obtain a reliable and powerful unit in a solid configuration at a very attractive price. This heating oven is preferred for tasks in basic research, as well as in human and veterinary medicine.





E 28 model

BENEFITS

- Uniform drying conditions thanks to APT.line[™] technology
- Identical test conditions throughout the chamber interior independent of sample size and quantity
- Outstanding thermal insulation saves operation costs

MAIN FEATURES

- Temperature range: +60 °C to +230 °C
- Adjustable exhaust air flap
- Hydro-mechanical thermostat
- Class 1 temperature limitation device (Model E028-230V-T)
- Timer 0 120 min
- 2 chrome-plated racks

AVAILABLE SIZES

Model	Interior volume [L]	Online data sheets	
E 28	28	› go2binder.com/en-E28-CL	_

G

OPTIONAL EQUIPMENT

Access port with silicone plug	For introducing external measuring instruments. Large selection of different diameters and positions (left, right, or top)	
Viewing window and interior lighting	Viewing window in the door with interior lighting for easy samples viewing.	-
Shelf racks	Perforated shelves or shelf racks, chrome-plated or made from stainless steel.	
Services	You can access our installation, maintenance, calibration, validation, or warranty services by taking out a Maintenance Service Agreement or by simply getting in touch when you need us. More information is available in the accessories and services chapter.	
BINDER INDIVIDUAL	Using a BINDER series production product as a basis, we develop individual solutions as a single unit or in small batches. All products are tested, certified, and supplied with a full functional guarantee and warranty. More information is available in the accessories and services chapter.	-
► All extras online go2bind	er.com/en-options	


Description	E 28
Measures	
Interior volume [L]	28
Net weight of the unit (empty) [kg]	22
Permitted load [kg]	25
Load per rack [kg]	10
Wall clearance back [mm]	160
Wall clearance sidewise [mm]	100
Internal Dimensions	
Width [mm]	400
Height [mm]	280
Depth [mm]	250
Housing dimensions not incl. fittings and connections	
Width net [mm]	580
Height net [mm]	405
Depth net [mm]	425
Performance Data Temperature	
Temperature range [°C]	60230
Heating-up time to 150 °C [min]	36
Recovery time after 30 seconds door open at 150 °C [min]	19
Electrical data	
Nominal power (version 120 V) [kW]	0.8
Nominal power (version 230 V) [kW]	0.8
Fixtures	
Number of shelves (std./max.)	2/4

ORDERING INFORMATION [Art. No.]

Further details: Go to www.binder-world.com > Search > Enter Article Number

Voltage	Option model	E 28
230 V 1~ 50/60 Hz	Standard	9010-0001
120 V 1~ 60 Hz	Standard	9010-0106

DIMENSIONS incl. fittings and connections [mm]

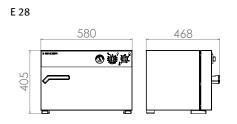

ORDERING INFORMATION OPTIONS

Designation	Description	E 28	Art. No.
Calibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test tempera- ture	•	8012-1114
Calibration certificate, temperature	for temperature, measurement in center of chamber at specified temperature	•	8012-1133

Designation	Description	E 28	Art. No.
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	•	8012-2250
Rack	chrome plated	•	8012-2032
Shelf, perforated	Stainless steel	•	8012-2160

Series ED Classic.Line | with natural convection

The strengths of a BINDER ED series drying chamber include routine drying and sterilization tasks up to +300 °C. Thanks to the natural convection, all thermal processes which use this drying chamber are highly efficient. The ED series ensures fast, uniform drying.





BD 400 model

BENEFITS

- Uniform drying conditions thanks to APT.line™ technology
- Identical test conditions throughout the chamber interior independent of sample size and quantity
- Outstanding thermal insulation saves operation costs

MAIN FEATURES

- Temperature range: +5 °C above ambient temperature to +300 °C
- APT.line[™] preheating chamber technology
- Natural convection
- Adjustable exhaust air flap
- Controller with timer functions
- 2 chrome-plated racks
- Class 2 independent adjustable temperature safety device (DIN 12880) with visual alarm

AVAILABLE SIZES

Model	Interior volume [L]	Online data sheets
ED 23	20	› go2binder.com/en-ED23-CL
ED 400	400	> go2binder.com/en-ED400-CL

G

OPTIONAL EQUIPMENT


Access port with silicone plug	For introducing external measuring instruments. Large selection of different diameters and positions (left, right, or top)	
Viewing window and interior lighting	Viewing window in the door with interior lighting for easy samples viewing.	
Shelf racks	Perforated shelves or shelf racks, chrome-plated or made from stainless steel.	
Services	You can access our installation, maintenance, calibration, validation, or warranty services by taking out a Maintenance Service Agreement or by simply getting in touch when you need us. More information is available in the accessories and services chapter.	
BINDER INDIVIDUAL	Using a BINDER series production product as a basis, we develop individual solutions as a single unit or in small batches. All products are tested, certified, and supplied with a full functional guarantee and warranty. More information is available in the accessories and services chapter.	-
► All extras online go2bind	er.com/en-options	

Description	ED 23	ED 400
Measures		
Interior volume [L]	20	400
Net weight of the unit (empty) [kg]	27	125
Permitted load [kg]	25	90
Load per rack [kg]	12	35
Wall clearance back [mm]	100	160
Wall clearance sidewise [mm]	100	100
Internal Dimensions		
Width [mm]	222	1000
Height [mm]	330	800
Depth [mm]	300	520
Housing dimensions not incl. fittings and connections		
Width net [mm]	435	1235
Height net [mm]	495	1025
Depth net [mm]	520	765
Performance Data Temperature		
Temperature range +5 °C above ambient temperature to [°C]	300	300
Temperature variation at 150 °C [± K]	2.8	3
Temperature fluctuation at 150 °C [± K]	0.5	0.5
Heating-up time to 150 °C [min]	28	71
Recovery time after 30 seconds door open at 150 °C [min]	5	31
Electrical data		
Nominal power (version 120 V) [kW]	0.8	-
Nominal power (version 208 V) [kW]	-	3.4
Nominal power (version 230 V) [kW]	0.8	-
Nominal power (version 400 V) [kW]	-	3.4
Fixtures		
Number of shelves (std./max.)	2/4	2/10

ORDERING INFORMATION [Art. No.] Further details: Go to www.binder-world.com > Search > Enter Article Number

Voltage	Option model	ED 23	ED 400
400 V 3~ 50/60 Hz	with RS 422	_	9010-0075
230 V 1~ 50/60 Hz	Standard	9010-0190	-
208 V 3~ 60 Hz	with RS 422	-	9010-0168
120 V 1~ 60 Hz	Standard	9010-0192	-

DIMENSIONS incl. fittings and connections [mm]

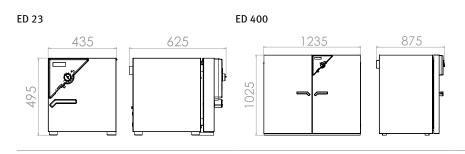

ORDERING INFORMATION OPTIONS

Designation	Description	ED 23	ED 400	Art. No.
	back			
	10 mm	•	•	8012-1282
	30 mm	•	•	8012-0540
	50 mm	•	•	8012-054
	left			
	10 mm	•	•	8012-127
	30 mm	•	•	8012-053
	50 mm	•	•	8012-054
	100 mm	-	•	8012-072
Access port with sili- cone plug	right			
	10 mm	•	•	8012-126
	30 mm	•	•	8012-053
	50 mm	•	•	8012-054
	100 mm	-	•	8012-072
	top			
	10 mm	•	•	8012-126
	30 mm	•	•	8012-053
	50 mm	•	•	8012-054
	100 mm	•	•	8012-054
Alarm function for	Switchable acoustic alarm, with adjustable limit value on the independent temperature safety device			
overheating	120/208 V option model	•	•	8012-130
	230/400 V option model	•	•	8012-129
Analog output 4-20 mA	for temperature values (output not adjustable)	•	•	8012-043
Calibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test temperature	•	•	8012-111
	for temperature, measurement in center of chamber at specified temperature	•	•	8012-113
Calibration certificate,	temperature measurement incl. certificate and 27 measuring points at specified tempera- ture	•	•	8012-159
temperature	temperature measurement incl. certificate, 15-18 measuring points at specified tempera- ture	•	•	8012-157
	temperature measurement incl. certificate, 9 measuring points at specified temperature	•	•	8012-155
Class 3.1 independent temperature safety device	with visual alarm (DIN 12880)	•	•	8012-003
		•	_	8012-049
)oor gasket	made of FKM, silicone-free	_	•	8012-049

Designation	Description	ED 23	ED 400	Art. No.
Door lock	lockable door handle	•	-	8012-1246
		-	•	8012-1253
Inner chamber, rein- forced	including 2 reinforced racks	-	•	8012-173
Measurement of air exchange rate	in accordance with ASTM D5374, definition and protocol according to ambient temperature	•	•	8012-019
Pt 100 temperature sensor	additional flexible Pt 100, interior, with external connection incl. LEMO connector (3-pin)	•	•	8012-117
	viewing window in door, 180 x 180 mm, and 15 W interior lighting			
	120 V option model	•	-	8012-037
Viewing window and	230 V option model	•	-	8012-037
interior lighting	viewing window in each door, 290 x 470 mm, and 30 W interior lighting			
	208 V option model	-	•	8012-027
	400 V option model	_	•	8012-026

Designation	Description	ED 23	ED 400	Art. No.
АРТ-СОМ™ 4 ВА-	for simple logging and documentation requirements with up to 5 networked units.			
SIC-Edition	version 4, BASIC edition	-	•	9053-0039
APT-COM™ 4 GLP-Editi-	for working under GLP-compliant conditions. Measured values are documented in a tam- per-proof way in line with the requirements of FDA Regulation 21 CFR 11.			
on	version 4, GLP edition	-	•	9053-004
APT-COM™ 4 PROFES-	convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units.			
SIONAL-Edition	version 4, PROFESSIONAL edition	-	•	9053-004
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	•	•	8012-225
	IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature, CO_2 , O_2 , pressure, depending on unit			
	Digital in PDF format	•	•	7057-000
Qualification docu-	Hard copy inside folder	•	•	7007-000
ments	IQ/OQ/PQ documents – supporting documents for validation performed by customers, according to customer requirements, PQ section added to qualification folder IQ/OQ; parameters: temperature, CO_2 , O_2 – or pressure, depending on unit			
	Digital in PDF format	•	•	7057-000
	Hard copy inside folder	•	•	7007-000
		•	_	8012-203
	chrome plated	-	•	8012-203
Rack		•	-	8012-204
	stainless steel	-	•	8012-215
Rubber pads	set anti-slip feet	•	-	8012-203
Shelf, perforated	Stainless steel	-	•	8012-216
Table on castors	stable cart, casters with locking brakes, dimensions: W 1,300 x D 800 x H 780 mm	_	•	9051-001

Series FD Classic.Line | with forced convection

A BINDER FD series heating oven is always used when fast drying and sterilization is required. Thanks to its fully homogeneous temperature distribution, quick dynamics and powerful fan, this heating oven saves valuable time.



FD 23 model

BENEFITS

- Uniform drying conditions thanks to APT.line™ technology
- Identical test conditions throughout the chamber interior independent of sample size and quantity
- Outstanding thermal insulation saves operation costs



MAIN FEATURES

- Temperature range: +5 °C above ambient temperature to +300 °C
- APT.line™ preheating chamber technology
- Forced convection
- Adjustable exhaust air flap
- Controller with timer functions
- 2 chrome-plated racks
- Class 2 independent adjustable temperature safety device (DIN 12880) with visual alarm

AVAILABLE SIZES

Model	Interior volume [L]	Online data sheets	
FD 23	20	> go2binder.com/en-FD23-CL	



OPTIONAL EQUIPMENT

Access port with silicone plug	For introducing external measuring instruments. Large selection of different diameters and positions (left, right, or top)	
Viewing window and interior lighting	Viewing window in the door with interior lighting for easy samples viewing.	
Shelf racks	Perforated shelves or shelf racks, chrome-plated or made from stainless steel.	
Services	You can access our installation, maintenance, calibration, validation, or warranty services by taking out a Maintenance Service Agreement or by simply getting in touch when you need us. More information is available in the accessories and services chapter.	
BINDER INDIVIDUAL	Using a BINDER series production product as a basis, we develop individual solutions as a single unit or in small batches. All products are tested, certified, and supplied with a full functional guarantee and warranty. More information is available in the accessories and services chapter.	-
► All extras online go2bind	er.com/en-options	

Description	FD 23	
Measures		
Interior volume [L]	20	
Net weight of the unit (empty) [kg]	27	
Permitted load [kg]	25	
Load per rack [kg]	12	
Wall clearance back [mm]	160	
Wall clearance sidewise [mm]	100	
Internal Dimensions		
Width [mm]	222	
Height [mm]	330	
Depth [mm]	300	
Housing dimensions not incl. fittings and connections		
Width net [mm]	435	
Height net [mm]	495	
Depth net [mm]	520	
Performance Data Temperature		
Temperature range +5 °C above ambient temperature to [°C]	300	
Temperature variation at 150 °C [± K]	2.5	
Temperature fluctuation at 150 °C [± K]	0.3	
Heating-up time to 150 °C [min]	25	
Recovery time after 30 seconds door open at 150 °C [min]	6	
Electrical data		
Nominal power (version 120 V) [kW]	0.8	
Nominal power (version 230 V) [kW]	0.8	
Fixtures		
Number of shelves (std./max.)	2/4	


ORDERING INFORMATION [Art. No.]

Further details: Go to www.binder-world.com > Search > Enter Article Number

Voltage	Option model	FD 23
230 V 1~ 50/60 Hz	Standard	9010-0194
120 V 1~ 60 Hz	Standard	9010-0196

DIMENSIONS incl. fittings and connections [mm]



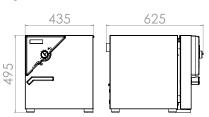

ORDERING INFORMATION OPTIONS

Designation	Description	FD 23	Art. No.
	back		
	10 mm	•	8012-128
	left		
	10 mm	•	8012-127
	30 mm	•	8012-135
	50 mm	•	8012-137
Access port with sili-	right		
cone plug	10 mm	•	8012-126
	30 mm	•	8012-135
	50 mm	•	8012-137
	top		
	10 mm	•	8012-126
	30 mm	•	8012-136
	50 mm	•	8012-138
	Switchable acoustic alarm, with adjustable limit value on the independent temperature safety device		
Alarm function for	120 V option model	•	8012-130
overheating	230 V option model	•	8012-130
Analog output 4-20 mA	for temperature values (output not adjustable)	•	8012-043
Calibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test tempera- ture	•	8012-111
	for temperature, measurement in center of chamber at specified temperature	•	8012-113
Calibration certificate,	temperature measurement incl. certificate and 27 measuring points at specified temperature	•	8012-159
emperature	temperature measurement incl. certificate, 15- 18 measuring points at specified temperature	•	8012-157
	temperature measurement incl. certificate, 9 measuring points at specified temperature	•	8012-155
Class 3.1 independent temperature safety device	with visual alarm (DIN 12880)	•	8012-114
Door gasket	made of FKM, silicone-free	•	8012-130
Door lock	lockable door handle	•	8012-124
Measurement of air exchange rate	in accordance with ASTM D5374, definition and protocol according to ambient temperature	•	8012-119
Pt 100 temperature sensor	additional flexible Pt 100, interior, with external connection incl. LEMO connector (3-pin)	•	8012-117
	viewing window in door, 180 x 180 mm, and 15 W interior lighting		
Viewing window and Interior lighting	120 V option model	•	8012-126
interior ugntilig	230 V option model	•	8012-125

Designation	Description	FD 23	Art. No.
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	•	8012-2250

Designation	Description	FD 23	Art. No.
	IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature, CO_2 , O_2 , pressure, depending on unit		
	Digital in PDF format	•	7057-0001
Qualification docu-	Hard copy inside folder	•	7007-0001
ments	IQ/OQ/PQ documents – supporting documents for validation performed by customers, according to customer requirements, PQ section added to qualification folder IQ/OQ; parameters: temperature, CO_2 , O_2 – or pressure, depending on unit		
	Digital in PDF format	•	7057-0005
	Hard copy inside folder	•	7007-0005
	chrome plated	•	8012-2031
Rack	stainless steel	•	8012-2046
Rack accessories	fasteners (1 set of 4) for additional security of racks	•	8012-0531
Shelf, perforated	Stainless steel	•	8012-2164

Series FED Classic.Line | with forced convection and enhanced timer functions

BINDER heating chambers of the FED series provide almost unlimited capacity and are particularly adaptable. With advanced timing functions and a adjustable fan, temperature and convection conditions are easily controlled.





BF 400 model

BENEFITS

- Uniform drying conditions thanks to APT.line[™] technology
- Identical test conditions throughout the chamber interior independent of sample size and quantity
- Outstanding thermal insulation saves operation costs

MAIN FEATURES

- Temperature range: +5 °C above ambient temperature to +300 °C
- APT.line[™] preheating chamber technology
- Adjustable fan speed
- Adjustable exhaust air flap
- Controller with expanded timer functions
- 2 chrome-plated racks
- Class 2 independent adjustable temperature safety device (DIN 12880) with visual alarm
- Computer interface: RS 422

AVAILABLE SIZES

Model	Interior volume [L]	Online data sheets
FED 400	400	› go2binder.com/en-FED400-CL

G

OPTIONAL EQUIPMENT

Access port with silicone plug	For introducing external measuring instruments. Large selection of different diameters and positions (left, right, or top)	
Viewing window and interior lighting	Viewing window in the door with interior lighting for easy samples viewing.	
Shelf racks	Perforated shelves or shelf racks, chrome-plated or made from stainless steel.	
Services	You can access our installation, maintenance, calibration, validation, or warranty services by taking out a Maintenance Service Agreement or by simply getting in touch when you need us. More information is available in the accessories and services chapter.	
BINDER INDIVIDUAL	Using a BINDER series production product as a basis, we develop individual solutions as a single unit or in small batches. All products are tested, certified, and supplied with a full functional guarantee and warranty. More information is available in the accessories and services chapter.	-
► All extras online go2bind	er.com/en-options	

Description	FED 400	
Measures		
Interior volume [L]	400	
Net weight of the unit (empty) [kg]	145	
Permitted load [kg]	90	
Load per rack [kg]	35	
Wall clearance back [mm]	160	
Wall clearance sidewise [mm]	100	
Internal Dimensions		
Width [mm]	1000	
Height [mm]	800	
Depth [mm]	510	
Housing dimensions not incl. fittings and connections		
Width net [mm]	1235	
Height net [mm]	1025	
Depth net [mm]	765	
Performance Data Temperature		
Temperature range +5 °C above ambient temperature to [°C]	300	
Temperature variation at 150 °C [± K]	3.8	
Temperature fluctuation at 150 °C [± K]	0.7	
Heating-up time to 150 °C [min]	29	
Recovery time after 30 seconds door open at 150 °C [min]	6	
Electrical data		
Nominal power (version 208 V) [kW]	3.4	
Nominal power (version 400 V) [kW]	3.4	
Fixtures		
Number of shelves (std./max.)	2/10	


ORDERING INFORMATION [Art. No.]

Further details: Go to www.binder-world.com > Search > Enter Article Number

Voltage	Option model	FED 400
400 V 3~ 50/60 Hz	Standard	9010-0216
208 V 3~ 60 Hz	Standard	9010-0217

DIMENSIONS incl. fittings and connections [mm]

FED 400

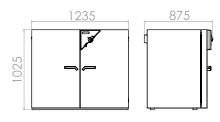

ORDERING INFORMATION OPTIONS

Designation	Description	FED 400	Art. No.
	back		
	10 mm	•	8012-1285
	30 mm	•	8012-1369
	50 mm	•	8012-1385
	left		
	10 mm	•	8012-1278
	30 mm	•	8012-136
	50 mm	•	8012-137
	100 mm	•	8012-150
Access port with sili-	right		
one plug	10 mm	•	8012-127
	30 mm	•	8012-135
	50 mm	•	8012-137
	100 mm	•	8012-150
	top		
	10 mm	•	8012-126
	30 mm	•	8012-136
	50 mm	•	8012-138
	100 mm	•	8012-138
			8012-138
larm function for	Switchable acoustic alarm, with adjustable limit value on the independent temperature safety device	•	0012 120
overheating	208 V option model	•	8012-130
	400 V option model	•	8012-130
Analog output 4-20 mA	for temperature values (output not adjustable)	•	8012-043
alibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test tempera- ture	•	8012-111
	for temperature, measurement in center of chamber at specified temperature	•	8012-113
Calibration certificate,	temperature measurement incl. certificate and 27 measuring points at specified temperature	•	8012-159
emperature	temperature measurement incl. certificate, 15-18 measuring points at specified temperature	•	8012-157
	temperature measurement incl. certificate, 9 measuring points at specified temperature	•	8012-155
Class 3.1 independent emperature safety levice	with visual alarm (DIN 12880)	•	8012-114
Door gasket	made of FKM, silicone-free	•	8012-131
)oor lock	lockable door handle	•	8012-125
	to increase air change rates		
an, reinforced	400 V option model	•	8012-020
IEPA air filter	On fresh air supply for unit; filter classification H14 (acc. to. EN 1822-1:2009, > 99.995% at 0.3 µm); not in conjunction with reinforced heating output, reinforced fan	•	8012-007
nner chamber, rein- orced	including 2 reinforced racks	•	8012-033
Measurement of air exchange rate	in accordance with ASTM D5374, definition and protocol according to ambient temperature	•	8012-119
Pt 100 temperature sensor	additional flexible Pt 100, interior, with external connection incl. LEMO connector (3-pin)	•	8012-117
	viewing window in each door, 290 x 470 mm, and 30 W interior lighting		
/iewing window and	208 V option model	•	8012-124
nterior lighting	400 V option model	•	8012-123

Designation	Description	FED 400	Art. No.
АРТ-СОМ™ 4 ВА-	for simple logging and documentation requirements with up to 5 networked units.		
SIC-Edition	version 4, BASIC edition	•	9053-0039
APT-COM™ 4 GLP-Editi-	for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11.		
on	version 4, GLP edition	•	9053-004
APT-COM™ 4 PROFES- SIONAL-Edition	convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units.		
SIONAL-Edition	version 4, PROFESSIONAL edition	•	9053-004
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	•	8012-2250
	IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature, CO_2 , O_2 , pressure, depending on unit		
	Digital in PDF format	•	7057-000
Qualification docu-	Hard copy inside folder	•	7007-000
ments	$IQ/OQ/PQ$ documents – supporting documents for validation performed by customers, according to customer requirements, PQ section added to qualification folder IQ/OQ; parameters: temperature, CO_2 , O_2 – or pressure, depending on unit		
	Digital in PDF format	•	7057-000
	Hard copy inside folder	•	7007-000
	chrome plated	•	8012-2036
Rack	stainless steel	•	8012-2159
Rack accessories	fasteners (1 set of 4) for additional security of racks	•	8012-053
Rack, reinforced	stainless steel, with fasteners (1 set of 4)	•	8012-034
Shelf, perforated	Stainless steel	•	8012-2162
Table on castors	stable cart, casters with locking brakes, dimensions: W 1,300 x D 800 x H 780 mm	•	9051-001

Series M Classic.Line | with forced convection and advanced program

functions

With a maximum of 300 °C and various programming options the M series material test chambers are ideal for material and aging tests. With a high air flow, the powerful fans ensure rapid heating.



M 53 model

BENEFITS

- Uniform drying conditions thanks to APT.line™ technology
- Identical test conditions throughout the chamber interior independent of sample size and quantity
- Outstanding thermal insulation saves operation costs
- Advanced programming options



MAIN FEATURES

- Temperature range: +5 °C above ambient temperature to +300 °C
- APT.line[™] preheating chamber technology
- Adjustable fan speed
- Program-controlled ventilation flap
- Color LCD controller with time-segment programming
- 2 chrome-plated racks
- Class 2 independent adjustable temperature safety device (DIN 12880) with visual alarm
- Computer interface: RS 422

AVAILABLE SIZES

Interior volume [L]	Online data sheets
53	› go2binder.com/en-M53-CL
115	> go2binder.com/en-M115-CL
240	> go2binder.com/en-M240-CL
400	> go2binder.com/en-M400-CL
720	> go2binder.com/en-M720-CL
	volume [L] 53 115 240 400




OPTIONAL EQUIPMENT


Access port with silicone plug	For introducing external measuring instruments. Large selection of different diameters and positions (left, right, or top)	
Viewing window and interior lighting	Viewing window in the door with interior lighting for easy samples viewing.	
Shelf racks	Perforated shelves or shelf racks, chrome-plated or made from stainless steel.	
Services	You can access our installation, maintenance, calibration, validation, or warranty services by taking out a Maintenance Service Agreement or by simply getting in touch when you need us. More information is available in the accessories and services chapter.	
BINDER INDIVIDUAL	Using a BINDER series production product as a basis, we develop individual solutions as a single unit or in small batches. All products are tested, certified, and supplied with a full functional guarantee and warranty. More information is available in the accessories and services chapter.	- 22
► All extras online go2bind	er.com/en-options	

Description	M 53	M 115	M 240	M 400	M 720
Measures					
Interior volume [L]	53	115	240	400	720
Net weight of the unit (empty) [kg]	61	89	131	173	203
Permitted load [kg]	40	50	70	90	120
Load per rack [kg]	15	20	30	35	45
Wall clearance back [mm]	160	160	160	160	160
Wall clearance sidewise [mm]	100	100	100	100	100
Internal Dimensions					
Width [mm]	400	600	800	1000	1000
Height [mm]	400	480	600	800	1200
Depth [mm]	340	410	510	510	610
Housing dimensions not incl. fittings and connections					
Width net [mm]	635	835	1035	1235	1235
Height net [mm]	780	865	985	1190	1690
Depth net [mm]	575	645	745	795	865
Performance Data Temperature					
Temperature range +5 °C above ambient temperature to [°C]	300	300	300	300	300
Temperature variation at 150 °C [± K]	1.3	2.5	1.5	1.5	1.9
Temperature fluctuation at 150 °C [± K]	0.3	0.3	0.3	0.3	0.3
Heating-up time to 150 °C [min]	15	16	19	18	21
Recovery time after 30 seconds door open at 150 °C [min]	3	3	3	3	3
Tabellen					
Air change (approx.) at 150 °C [x/h]	192	96	60	54	36
Electrical data					
Nominal power (version 230 V) [kW]	1.2	1.6	2.7	-	-
Nominal power (version 400 V) [kW]	-	-	-	3.4	5
Fixtures					
Number of shelves (std./max.)	2/5	2/6	2/7	2/10	2/16

ORDERING INFORMATION [Art. No.] Further details: Go to www.binder-world.com > Search > Enter Article Number


Voltage	Option model	M 53	M 115	M 240	M 400	M 720
400 V 3~ 50/60 Hz	Standard	-	-	-	9010-0204	9010-0205
230 V 1~ 50/60 Hz	Standard	9010-0201	9010-0202	9010-0203	-	-

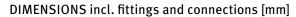

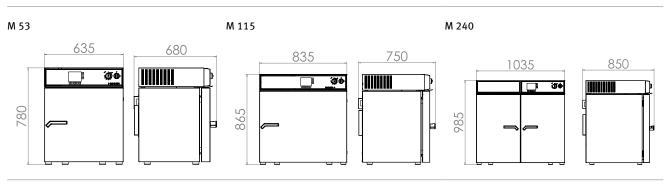

M 400

M 720

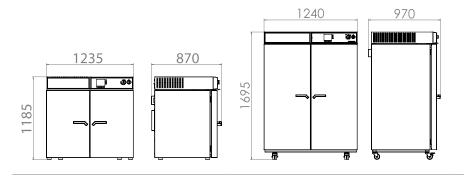

ORDERING INFORMATION OPTIONS

Designation	Description	M 53	M 115	M 240	M 400	M 720	Art. No.
	back						
	10 mm	•	•	•	•	•	8012-1287
	30 mm	•	•	•	•	•	8012-137
	50 mm	•	•	•	•	•	8012-138
	left						
	10 mm	•	•	٠	•	•	8012-128
	30 mm	•	•	•	•	•	8012-136
	50 mm	•	•	•	•	•	8012-137
Access port with sili-	100 mm	•	•	•	•	•	8012-151
cone plug	right						
	10 mm	•	•	•	•	•	8012-127
	30 mm	•	•	•	•	•	8012-135
	50 mm	•	•	•	•	•	8012-137
	100 mm	•	•	•	•	•	8012-150
	top						
	10 mm	•	•	•	•	•	8012-126
	30 mm	•	•	•	•	•	8012-136
	50 mm	•	•	•	•	•	8012-138
Analog output 4-20 mA	for temperature values (output not adjustable)	•	•	•	•	•	8012-048
Calibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test temperature	•	•	٠	•	•	8012-112

Designation	Description	M 53	M 115	M 240	M 400	M 720	Art. No.
	for temperature, measurement in center of chamber at speci- fied temperature	•	•	•	•	•	8012-1142
Calibration certificate,	temperature measurement incl. certificate and 27 measuring points at specified temperature	•	•	•	•	•	8012-160
emperature	temperature measurement incl. certificate, 15- 18 measuring points at specified temperature	•	•	•	•	•	8012-158
	temperature measurement incl. certificate, 9 measuring points at specified temperature	•	•	•	•	•	8012-155
		•	-	-	-	-	8012-130
		-	•	-	-	-	8012-130
Door gasket	made of FKM, silicone-free	-	-	•	-	-	8012-130
		-	-	-	•	-	8012-131
		-	-	_	-	•	8012-131
		•	•	_	-	_	8012-125
Door lock	lockable door handle	-	_	•	•	•	8012-125
Gas-tight option model	Incl. additional measures for the greatest possible impermea- bility to gas; not in conjunction with access port or interior lighting option	•	•	_	_	_	8012-018
IEPA air filter	On fresh air supply for unit; filter classification H14 (acc. to. EN 1822-1:2009, > 99.995% at 0.3 μm); not in conjunction with reinforced heating output, reinforced fan	•	•	•	•	•	8012-022
nert gas connection	with gas inlet and outlet, Ø 10 mm, incl. additional measures to provide the greatest possible impermeability to gas; not in conjunction with access port or interior lighting option	•	•	_	_	_	8012-154
	including 2 reinforced racks						
nner chamber, rein-		-	_	•	-	-	8012-034
nner chamber, rein- orced	max. total load 250 kg	-	-	-	•	-	8012-034
			_	-	-	•	8012-034
Lockable controller keyboard	Lockable controller keyboard with key-operated switch	•	•	•	•	•	8012-052
Measurement of air exchange rate	in accordance with ASTM D5374, definition and protocol according to ambient temperature	•	•	•	•	•	8012-119
Pt 100 temperature sensor	additional flexible Pt 100, interior, for displaying the tempera- ture on the unit display	•	•	•	•	•	8012-033
	viewing window in door, 230 x 230 mm, and 15 W interior lighting						
	230 V option model	•	-	-	-	-	8012-005
	viewing window in door, 320 x 260 mm, and 15 W interior lighting						
	230 V option model	-	•	-	-	-	8012-005
/iewing window and nterior lighting	viewing window in each door, 290 x 470 mm, and 30 W interior lighting						
		-	-	-	•	-	8012-006
	400 V option model	_	_	_	-	•	8012-006
	viewing window in each door, 350 x 240 mm, and 30 W interior lighting						
	230 V option model	_	_	•		_	8012-005

Designation	Description	M 53	M 115	M 240	M 400	M 720	Art. No.
APT-COM™ 4 BA-	for simple logging and documentation requirements with up to 5 networked units.						
SIC-Edition	version 4, BASIC edition	•	•	•	•	•	9053-003
APT-COM™ 4 GLP-Editi- on	for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11.						
	version 4, GLP edition	•	•	•	•	•	9053-004
APT-COM™ 4 PROFES- SIONAL-Edition	convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units.						
	version 4, PROFESSIONAL edition	•	•	•	•	•	9053-004
Base	the base equipped with casters is used for the secure positio-	-	•	-	-	-	9051-003
	ning and leveling of a BINDER drying and heating chamber	-	-	•	-	-	9051-003
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	•	•	•	•	•	8012-225
	IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature, CO ₂ , O ₂ , pressure, depending on unit						
Qualification docu- ments	Digital in PDF format	•	•	•	•	•	7057-000
	Hard copy inside folder	•	•	•	•	•	7007-000
	IQ/OQ/PQ documents – supporting documents for validation performed by customers, according to customer requirements, PQ section added to qualification folder IQ/OQ; parameters: temperature, CO ₂ , O ₂ – or pressure, depending on unit						
	Digital in PDF format	•	•	•	•	•	7057-000
	Hard copy inside folder	•	•	•	•	•	7007-000
		•	-	-	-	-	8012-203
		-	•	-	-	-	8012-203
	chrome plated	-	-	•	-	-	8012-203
		-	-	-	•	-	8012-203
Rack		-	-	-	-	•	8012-203
nuen		•	-	-	-	-	8012-204
		-	•	-	-	-	8012-204
	stainless steel	-	-	•	-	-	8012-208
		-	-	-	•	-	8012-215
		-	-	-	-	•	8012-208
Rack accessories	fasteners (1 set of 4) for additional security of racks	•	•	•	•	•	8012-053
		-	-	•	-	-	8012-034
Rack, reinforced	stainless steel, with fasteners (1 set of 4)	-	-	-	•	-	8012-034
		-	-	-	-	•	8012-037
Rubber pads	set anti-slip feet	•	•	•	•	-	8012-203
		•	-	-	-	-	8012-208
		-	•	-	-	-	8012-216
Shelf, perforated	Stainless steel	-	-	•	-	-	8012-208
		-	-	-	•	-	8012-216
		-	-	-	-	•	8012-208
Table on castors	stable cart, casters with locking brakes, dimensions: W 1,000 x D 800 x H 780 mm	•	•	-	-	-	9051-001
	stable cart, casters with locking brakes, dimensions: W 1,300			•	•		9051-001

Drying and heating chambers Solid.Line

Series ED-S Solid.Line

Series FD-S Solid.Line

Page 66

Page 70

With its new Solid.Line range, BINDER is expanding its product portfolio of drying and heating chambers, incubators, and constant climate chambers. This new line embodies solid product design with great value for money. The new Solid.Line is the ideal addition to BINDER's popular premium products, thus complementing the entire product range with BINDER's quality guarantee – Made in Germany.

PRODUCT SELECTION DRYING AND HEATING CHAMBERS SOLID.LINE

Features	Series ED-S Solid.Line > Page 66	Series FD-S Solid.Line > Page 70
Model		
Model 56	Ð	٠
Model 115	•	•
General data		
Temperature range +7 °C above ambient tempe- rature to [°C]	250	_
Temperature range +10 °C above ambient tempe- rature to [°C]	_	250
•Standard equipment		

Standard equipment


ADVANTAGES AT A GLANCE



EASY OPERATION OF THE UNIT

Series ED-S, FD-S

Using controller with digital display in °C and °F.

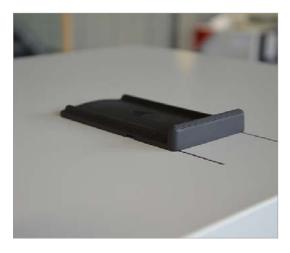

ADJUSTABLE EXHAUST AIR FLAP ABOVE THE UNIT

Series ED-S, FD-S For setting the air change rates.



COMPLETELY HOMOGENEOUS TEMPERATURE DISTRIBUTION

Series FD-S

This feature provides excellent spatial and temporal temperature accuracy with the all-new generation of APT.line™ technology.

Series ED-S Solid.Line | with natural convection

The new products in the Solid.Line impress with their tried-and-tested BINDER quality and reliability. They are suited to drying and heating. The attributes of the Solid.Line are geared towards the application in question and therefore used in research and quality assurance.



ED-S 56 model

BENEFITS

- Easy handling
- Accurate temperature control
- Simple and ergonomic door opening

MAIN FEATURES

- Temperature range: +7 °C above ambient temperature to +250 °C
- APT.line[™] preheating chamber technology
- Natural convection
- Adjustable exhaust air flap
- Controller with timer function
- 1 chrome-plated rack, incl. shelf supports
- Class 2 integrated independent adjustable temperature safety device (DIN 12880) with visual alarm

AVAILABLE SIZES

Model	Interior volume [L]	Online data sheets
ED-S 56	62	› go2binder.com/en-ED-S56
ED-S 115	118	› go2binder.com/en-ED-S115

G

OPTIONAL EQUIPMENT

Services

You can access our installation, maintenance, calibration, validation, or warranty services by taking out a Maintenance Service Agreement or by simply getting in touch when you need us. More information is available in the accessories and services chapter.



► All extras online go2binder.com/en-options


Description	ED-S 56	ED-S 115
Measures		
Interior volume [L]	62	118
Net weight of the unit (empty) [kg]	36	50
Permitted load [kg]	45	75
Load per rack [kg]	15	20
Wall clearance back [mm]	160	160
Wall clearance sidewise [mm]	100	100
Internal Dimensions		
Width [mm]	400	550
Height [mm]	440	550
Depth [mm]	350	390
Housing dimensions not incl. fittings and connections		
Width net [mm]	614	764
Height net [mm]	625	735
Depth net [mm]	577	616
Performance Data Temperature		
Temperature range +7 °C above ambient temperature to [°C]	250	250
Temperature variation at 150 °C [± K]	3.4	2.8
Temperature fluctuation at 150 °C [± K]	0.8	0.9
Heating-up time to 150 °C [min]	60	55
Electrical data		
Nominal power (version 120V) [kW]	1.15	1.35
Nominal power (version 230V) [kW]	1.05	1.25
Environment-specific data		
Energy consumption at 150 °C [Wh/h]	210	300
Fixtures		
Number of shelves (std./max.)	1/3	1/5

ORDERING INFORMATION [Art. No.]

Further details: Go to www.binder-world.com > Search > Enter Article Number

Voltage	Option model	ED-S 56	ED-S 115
230 V 1~ 50/60 Hz	Standard	9090-0014	9090-0020
120 V 1~ 60 Hz	Standard	9090-0015	9090-0021

DIMENSIONS incl. fittings and connections [mm]

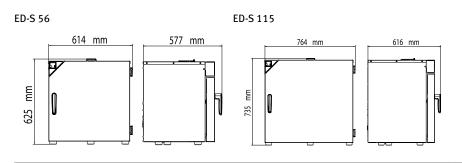

Description	ED-S 56	ED-S 115	Art. No.
Incl. shelf supports			
abrama alatad	•	-	8012-1807
chrome plated	-	•	8012-1808
concentrated, for gentle remove of residual contaminants; 1 kg	٠	•	8012-2250
set anti-slip feet	•	•	8012-2030
	Incl. shelf supports chrome plated concentrated, for gentle remove of residual contaminants; 1 kg	Incl. shelf supports • chrome plated • concentrated, for gentle remove of residual contaminants; 1 kg •	incl. shelf supports - chrome plated - concentrated, for gentle remove of residual contaminants; 1 kg •

Series FD-S Solid.Line | with forced convection

The new products in the Solid.Line impress with their tried-and-tested BINDER quality and reliability. They are suited to drying and heating. The attributes of the Solid.Line are geared towards the application in question and therefore used in research and quality assurance.



FD-S 115 model

BENEFITS

- Easy handling
- Accurate temperature control
- Simple and ergonomic door opening

MAIN FEATURES

- Temperature range: +10 °C above ambient temperature to +250 °C
- APT.line[™] preheating chamber technology
- Forced convection
- Adjustable exhaust air flap
- Controller with timer function
- 1 chrome-plated rack, incl. shelf supports
- Class 2 integrated independent adjustable temperature safety device (DIN 12880) with visual alarm

AVAILABLE SIZES

Model	Interior volume [L]	Online data sheets
FD-S 56	55	› go2binder.com/en-FD-S56
FD-S 115	106	› go2binder.com/en-FD-S115

G

OPTIONAL EQUIPMENT

Services

You can access our installation, maintenance, calibration, validation, or warranty services by taking out a Maintenance Service Agreement or by simply getting in touch when you need us. More information is available in the accessories and services chapter.



► All extras online go2binder.com/en-options


Description	FD-S 56	FD-S 115
Measures		
Interior volume [L]	55	106
Net weight of the unit (empty) [kg]	36	48
Permitted load [kg]	45	75
Load per rack [kg]	15	20
Wall clearance back [mm]	160	160
Wall clearance sidewise [mm]	100	100
Internal Dimensions		
Width [mm]	400	550
Height [mm]	440	550
Depth [mm]	310	350
Housing dimensions not incl. fittings and connections		
Width net [mm]	614	764
Height net [mm]	625	735
Depth net [mm]	577	616
Performance Data Temperature		
Temperature range +10 °C above ambient temperature to [°C]	250	250
Temperature variation at 150 °C [± K]	2.6	2.0
Temperature fluctuation at 150 °C [± K]	0.4	0.4
Heating-up time to 150 °C [min]	15	20
Tabellen		
Air change (approx.) at 100 °C [x/h]	4	2
Electrical data		
Nominal power (version 120V) [kW]	1.2	1.4
Nominal power (version 230V) [kW]	1.1	1.3
Environment-specific data		
Energy consumption at 150 °C [Wh/h]	290	340
Fixtures		
Number of shelves (std./max.)	1/3	1/5

ORDERING INFORMATION [Art. No.]

Further details: Go to www.binder-world.com > Search > Enter Article Number

Voltage	Option model	FD-S 56	FD-S 115
230 V 1~ 50/60 Hz	Standard	9090-0018	9090-0024
120 V 1~ 60 Hz	Standard	9090-0019	9090-0025

DIMENSIONS incl. fittings and connections [mm]

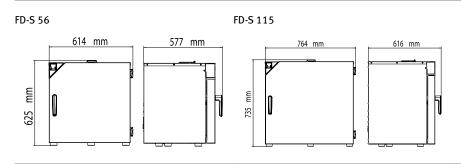

Designation	Description		FD-S	Art. No.
		56	115	
Grid shelf	Incl. shelf supports			
	chrome plated	•	-	8012-1807
		-	•	8012-1808
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg		•	8012-2250
Rubber pads	set anti-slip feet	٠	•	8012-2030

BINDER products

CO ₂ incubators				
Incubators				
Cooling incubators				
Growth chambers				
Ultralow temperature freezers				
Drying and heating chambers				
Safety drying chambers				
Vacuum drying chambers				
Constant climate chambers				
Dynamic climate chambers				
Battery test chambers				
Walk-in-chambers				
Multi Management Software				
Image: Second system Never miss a product novelty - Subscribe and keep up to speed with Subscribe and keep up to speed with Image: Second system Subscription				

go2binder.com/en-newsletter-subscription



Best conditions for your success

BINDER GmbH Tuttlingen, Germany

Tel +49 7462 2005 0 info@binder-world.com www.binder-world.com

BINDER Asia Pacific (Hong Kong) Ltd. Kowloon, Hong Kong, P.R. China

Tel +852 39070500 asia@binder-world.com

BINDER Environmental Testing Equipment (Shanghai) Co., Ltd. Shanghai, P.R.China

Tel +86 21 685 808 25 china@binder-world.com **BINDER Inc.** Bohemia, NY, USA

Tel +1 631 224 4340 usa@binder-world.com