



# Aeros<sup>®</sup>

Patent Pending

*Smart Color Measurement that Works for You*




measure  
color  
the way your  
eye sees it



# The *Smart* Way to Measure Color

Measuring smooth samples is a breeze for most color measurement instruments. But what about inconsistent, heavily textured, non-uniform samples like these?

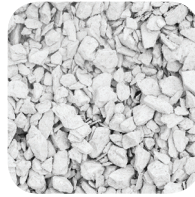

Rather than modifying your samples to work for your instrument, we built Aeros to work for your samples.

## Smart Sensor – Smart. Simple. Accurate.



### Auto-height Positioning

Positions the Smart Sensor to the ideal sample height and automatically adjusts for variations



### A 'One touch' out of the box color measurement

Reduced operator learning time, increases workflow and productivity

## Smart Sample platform - No Hassle. No Wasted Time.

Simply place the sample and hit measure



Samples automatically rotate beneath the Smart Sensor  
- 5 seconds, 35 measurements, 34.6 square inches of sample measured



Measures the largest sample area of any instrument for the most measurement accuracy and consistency available


## Smart Display – Durable Cognitive Computing

Designed for Industry 4.0

- Easy to use Hi-Res, Full Color Touch Interface
  - Durable industrial glass
  - Touch, Swipe, Type technology
  - Easy to read and view color data and graphics
- Embedded Color Quality Control Software
  - No PC required
  - All the color scales, indices and metrics you need
  - Intuitive and easy to learn user workflow



## Smart Communications – All in one technology

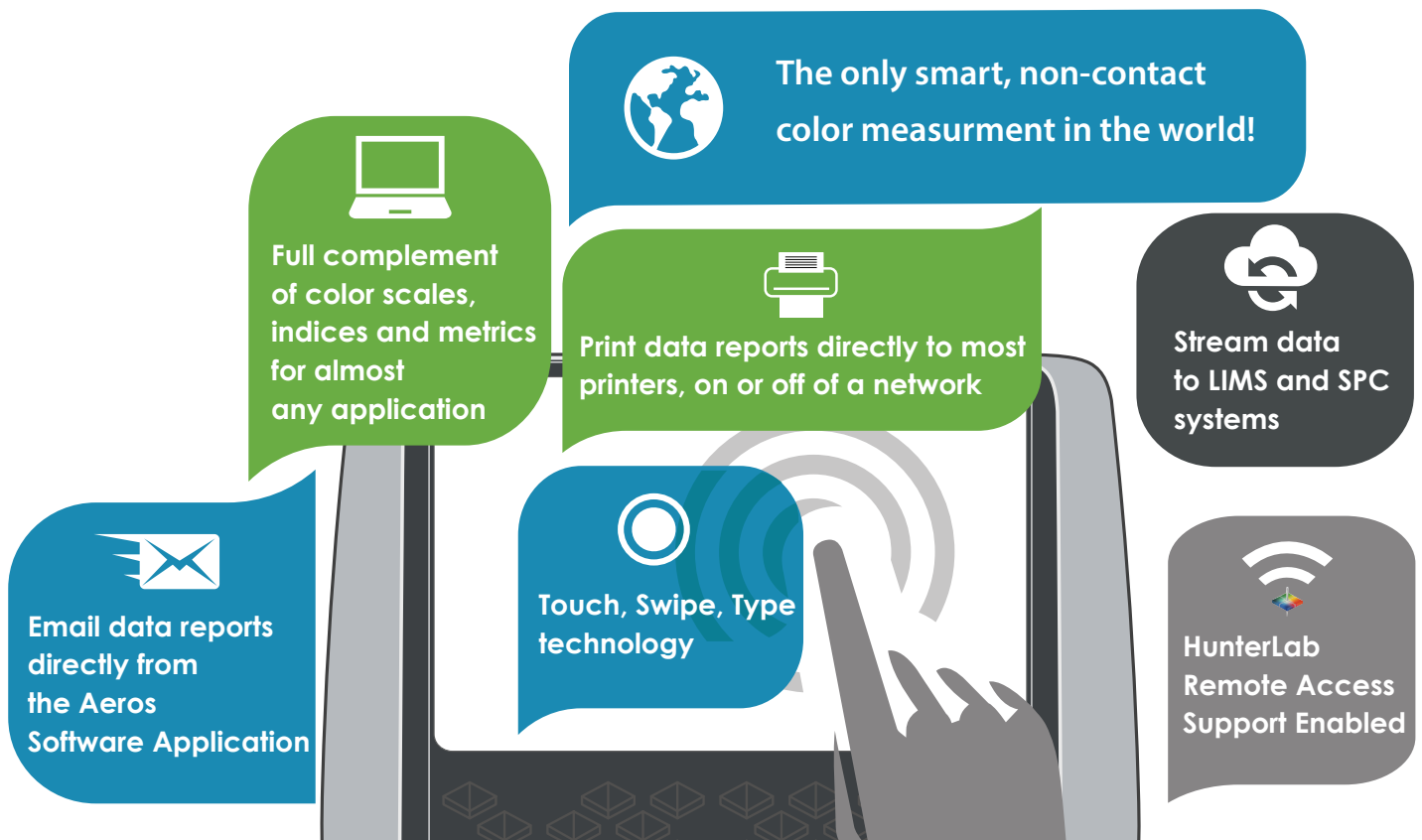



# Smart Color Measurement That Works for You

- Measure your samples as they are, in their natural state
- Fast, no worry sample preparation
- No need to measure through plastic or glass dishes
- No mess, minimal cleanup time
- Faster, easier, one button color measurement
- Accurate and consistent - true color results



*Designed for Industry 4.0*



## HunterLab Remote Access Support Enabled (Includes 1 year subscription)

HunterLab Remote Access Support is a unique support tool that enables the HunterLab support team to directly interface with your instrument to perform diagnostics, update instrument software and provide live operator training.

To learn more go to: <https://www.hunterlab.com/remote-access-support.html>



## What's in the box - *Everything you need to get it right.*



Aeros Reflectance Spectrophotometer



Standards box includes:

- Calibrated white tile (NIST Certificate of Traceability)
- Green diagnostic tile
- Black glass

Two Sample Trays:

- 30.5 cm (12 in) Dia x 2.54 cm (1 in) Depth
- 15.2 cm (6 in) Dia x 2.54 cm (1 in) Depth

Power cord (110/220 V)

For more information go to [www.hunterlab.com](http://www.hunterlab.com)  
or contact your local HunterLab representative.

HunterLab and Aeros (pat. pend.) are Trademarks of Hunter Associates Laboratory, Inc.  
Content subject to change without notice