

刀切式粉碎機 Cutting Mill



SW-2




CW-1



刀片內部
A close-up shot of the blade set



篩網組：
3 standard sieves included :
20 mesh x 1, 40 mesh x 1, 60 mesh x 1

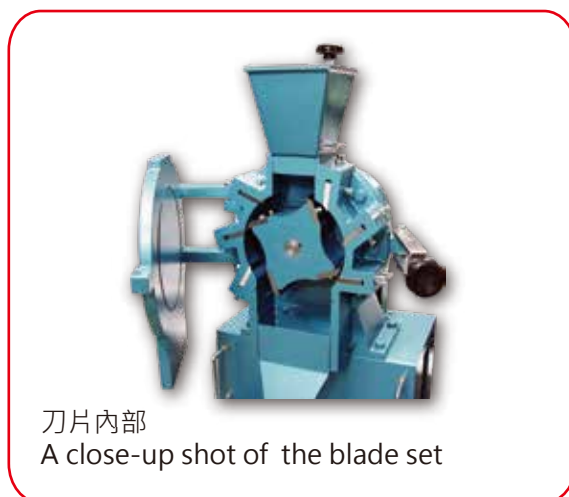

刀片內部
A close-up shot of the blade set

特點：

- 本機是採用鋼鐵製成，粉碎室週邊的固定刀及旋轉刀，將試料切斷粉碎且不減少水份，適合穀物、種子、纖維、皮革類之試料製作，粉碎刻度可由粉碎室下部之篩網控制。

Feature :

- Made of steel. The fixed knife and revolving Knife cut the materials into small pieces without reducing their water content. Applicable for grain, seeds, plastic fibers, and leathers. Size of the particles is controlled by the sieve under the mill pot.

Model	SW-2	CW-1
篩網孔徑 Sieve Opening	20, 40, 60 (mesh)	0.5, 1, 2 (mm)
入料口 Feeding Inlet Size	上 top \varnothing 6.2 cm / 下 bottom \varnothing 1.3 cm	160 x 160 mm
粉碎量 Capacity	5 kg / H	20 kg / H
馬力 Motor	200 W	750 W
尺寸 Size (W x L x H mm)	210 x 290 x 240	460 x 450 x 1100
重量 Weight	12.5 Kgs	127 Kgs
電源 Power	AC 110V / 220V, 50 / 60Hz	AC 220V, 1 Φ or 3 Φ